

Order ...

Institution / Company

Name / Position

Street

Zip Code / City

phone

fax

email

Date / Signature

Authorized h/p/cosmos dealer

h/p/cosmos sports & medical gmbh, Am Sportplatz 8, DE 83365 Nussdorf-Traunstein / Germany  
phone +49 / 8669 / 86 42 0, fax +49 / 8669 / 86 42 49, sales@h-p-cosmos.com, www.h-p-cosmos.com

number	order no.	item	price item	total price
_____	cos14825	h/p/cosmos sirius lactate test meter	_____	_____
_____	cos14854	h/p/cosmos sirius test strips, box 24 pcs.	_____	_____
_____	cos14827	h/p/cosmos sirius test strips, box 72 pcs.	_____	_____
_____	cos11350	Latex gloves, S, box 100 pieces	_____	_____
_____	cos11351	Latex gloves, M, box 100 pieces	_____	_____
_____	cos11355	Kidney dish, white	_____	_____
_____	cos11352	Desinfectant spray, 250 ml	_____	_____
_____	cos11652	Disposable lancets, sterile, 200 pieces	_____	_____
_____	cos11356	Plaster, Box 420 pieces	_____	_____
_____	cos11353	Swabs, 1000 pieces	_____	_____
_____	cos11930	Medi-Box (... for used lancets)	_____	_____
_____	cos14708	h/p/cosmos para analysis 2.0 pro	_____	_____
		price	_____	_____
		tax	_____	_____
		total price incl. tax	_____	_____

h/p/cosmos sirius®

lactate test meter

fast - precise - reliable



h/p/cosmos®

[cos15112] made in germany

© 06/2005 h/p/cosmos

id-card


## h/p/cosmos sirius®

The h/p/cosmos sirius lactate test meter always realises the demand for a fast and precise lactate analysis. With easy handling the menu navigation scrolling wheel supports the user in concentrating on the subject during the performance test.

The reliable measurement result appears online within 15 seconds.

Once the sufficient lactate value for the assesment of the subject has been reached, the test can be immediately ended without further strain on the subject.

As many as 250 results can be stored in the h/p/cosmos sirius and later transferred via interface into the desired test protocol of the software h/p/cosmos para analysis® 2.0 professional. This software offers a wide range of sports medical evaluations in performance diagnostics – from controlling ergometers to creating individual training programs.

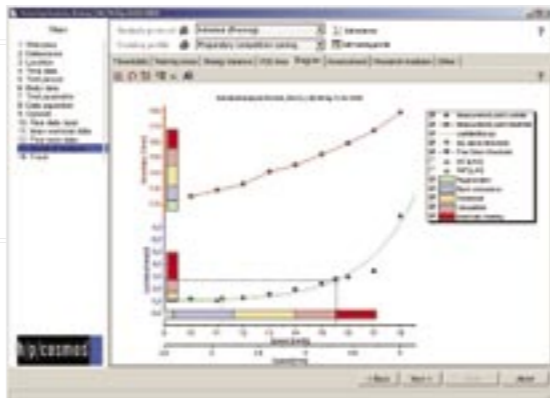

Measuring method:	enzymatical – amperometric
Sample quantity:	5 µl capillary blood
Measuring duration:	15 seconds
Measuring range:	0.5 ... 25 mmol/l
Measuring temperature:	min. 5° Celsius ... max. 45° Celsius
Humidity:	85 % relative humidity
Accuracy:	3 ... 8% depending on concentration
Memory capacity:	250 values with additional parameters
Weight:	approx. 80 g incl. batteries
Dimensions:	91 x 55 x 24 mm / 3.5 x 2.1 x 0.9 inches
Display:	LC-display with symbols
Power supply:	2 x 1.5 volt AAA/LR03 batteries
Interface:	serial interface / 2.5" cinch socket




Special test strip for precisely defined blood sample



Multi-functional scrolling wheel for easy handling

## h/p/cosmos sirius®

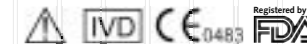

Manufacturer: SensLab GmbH  
Bautzner Straße 67, DE 04347 Leipzig



Interface for data transfer